

BLUES COMFORT

LINET

Comfort that adapts to every moment



Relaxing armchair and a comfortable extra bed in one.

Blues Comfort

It is ideal for maternity wards and children's wards, but also other facilities where a parent or partner stays overnight. It also serves to rest hard-working staff. It finds further use in nursing homes, where it helps people to share closeness and be with an elderly or ill person overnight at their bedside.

- The simple unfolding mechanism made of the metal construction ensures stability when lying on the armchair.
- The armchair can be upholstered in easy-to-clean fabric with special finishes, including antibacterial, UV-, abrasion- and chemical-resistant properties.
- The chair allows three relaxation positions - a sitting position, a relaxation position with the footrest pulled out and a sleeping position



Wooden visible elements including armrest in 9 shades of beech wood.

COLOUR CHART

Viny 1000	Viny 6275	Viny 7095	Viny 8067	Viny 8333
Viny 6056	Viny 2120	Viny 2216	Viny 5109	Viny 5203
Viny 4057	Viny 4697	Viny 1204	Viny 1999	

Blues Comfort Practic HPL



Solid plastic armrests are very durable and easy to maintain - for increased durability.



WOODEN DECORS

(only for version Blues Comfort)



ARCTIC	NATURE	HONEY	CATANO	TOBACCO
LARGO	WALNUT	BROWN	BLACK	



BASIC INFORMATION

- The strong, robust and stable folding mechanism ensures a comfortable and safe sleep. It consists of a galvanised metal structure which is protected by a high-quality powder-coated paint.
- We use welding technology to join the vinyl, which ensures a 100% hygienic and fluid-tight connection.
- The materials meet international standards for hygiene and fire resistance.
- BS 5852 Non-combustible
- DIN EN1021- 1 & 2
- Self-extinguishing
- Antibacterial
- Waterproof

TECHNICAL DATA

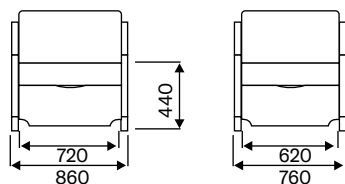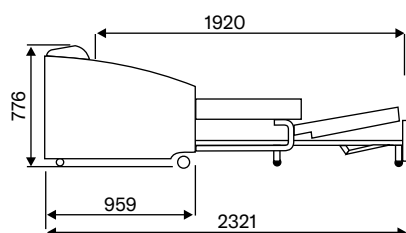

4 Wheels for better manipulation,
2 front wheels are 7,5cm with breaks
approximate weight 75kg
maximum weight of user 150kg



LINET

Members of LINET Group

LINET spol. s r.o.

Želečnice 5 | 274 01 Slaný | Czech Republic

tel.: +420 312 576 400 | fax: +420 312 522 668 | e-mail: info@linet.com | www.linet.com



www.linet.com